




Arrow®

NextStep® Antegrade Hemodialysis Catheter

The next step in dialysis catheter technology

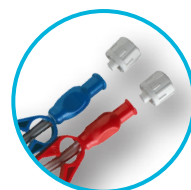

Arrow® NextStep® Hemodialysis Catheterization Kit

Step-Tip Ease Meets Split-Tip Flow

The Arrow® NextStep® Antegrade Catheter is a chronic hemodialysis catheter designed to combine:

- A step-tip catheter's ease of insertion
- A split-tip catheter's sustained high flow

Ease of Insertion

The Arrow® NextStep® Antegrade Catheter's unique tip is designed for:



- Smooth transition of over-the-wire insertion and exchanges
- Smooth transition through a sheath

Designed for sustained high flow and minimal recirculation

To reduce recirculation and deliver high flow, the Arrow® NextStep® Antegrade's tip has two unique, complementary features:

- The reversed ports are designed to take better advantage of blood flow dynamics
 - The venous port releases blood in the SVC
 - The arterial port draws blood from the right atrium
 - Side holes engineered to pull blood from all sides
- The ports are designed with significant separation to both enhance flow and minimize recirculation

Engineered For High Flow¹

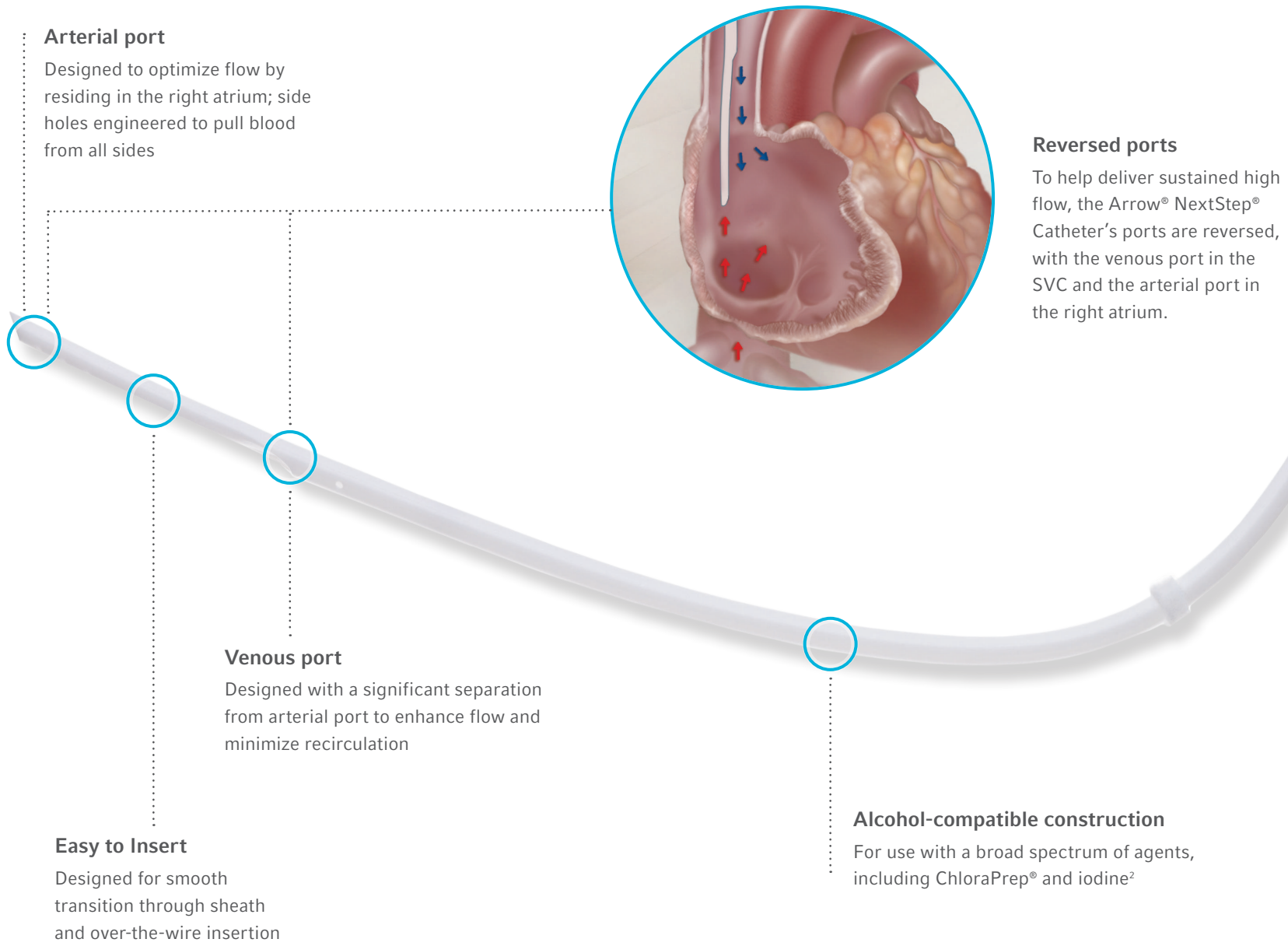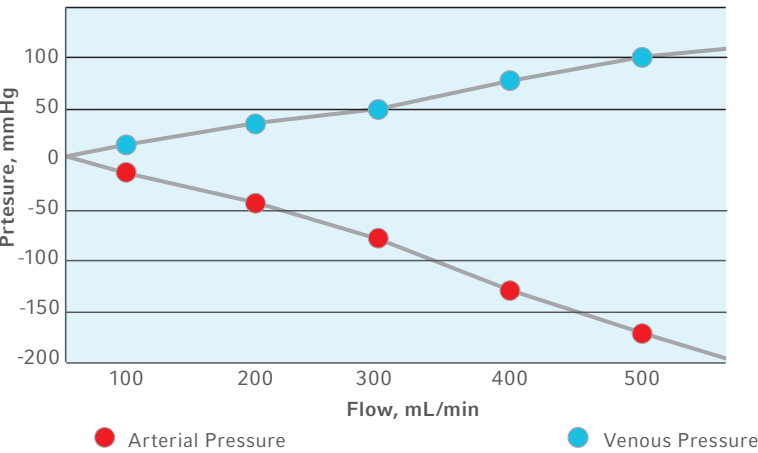

Arrow® NextStep® Hemodialysis Catheterization Kits

The Arrow® NextStep® Antegrade hemodialysis catheterization Kit includes:

- 1 Catheter: 15 Fr.
- 1 Arrow SmartSeal® Hemostatic Peelable Dialysis Sheath: 16 Fr.
- 1 Spring-Wire guide: .038" (.97 mm) x 39-1/2" (100 cm) straight stiff tip on one end—"J" tip on other with Arrow Advancer™
- 1 Introducer needle: 18 Ga. x 2-1/2" (6.35 cm)
- 1 Pre-loaded tunneler (metal) with Tunneling Sheath
- 2 Luer-Lock caps
- 1 Dressing: Tegaderm® 10 cm x 12 cm
- 1 Tissue dilator: 12 Fr.
- 1 Tissue dilator: 14 Fr.
- 1 SharpsAway II™ Locking Disposal Cup
- 1 Safety scalpel: #11
- 1 Patient ID card



Arrow® NextStep® Antegrade
Hemodialysis Catheterization Kit
(CS-15232-SFX)

To learn more about Arrow® NextStep® Catheters:
contact your Teleflex Representative

Ordering Information

Arrow® NextStep® Antegrade Chronic Hemodialysis Catheter Kits with SmartSeal® Hemostatic Sheath and Arrow Simplicity® Color-Coded Micro-Introducer

COLOR CODE	PRODUCT NO.	DESCRIPTION	TIP-TO-CUFF INSERTION LENGTH (CM)
19	CS-15192-SFXM	Arrow 15 Fr. x 7.5" (19 cm) Radiopaque Polyurethane Catheter with Color-Coded 21 Ga. Micro-Introducer Kit	19
23	CS-15232-SFXM	Arrow 15 Fr. x 9" (23 cm) Radiopaque Polyurethane Catheter with Color-Coded 21 Ga. Micro-Introducer Kit	23
27	CS-15272-SFXM	Arrow 15 Fr. x 10.5" (27 cm) Radiopaque Polyurethane Catheter with Color-Coded 21 Ga. Micro-Introducer Kit	27
31	CS-15312-SFXM	Arrow 15 Fr. x 12" (31 cm) Radiopaque Polyurethane Catheter with Color-Coded 21 Ga. Micro-Introducer Kit	31

Arrow® NextStep® Antegrade Chronic Hemodialysis Catheter Kits with SmartSeal® Hemostatic Dialysis Sheath

COLOR CODE	PRODUCT NO.	DESCRIPTION	TIP-TO-CUFF INSERTION LENGTH (CM)
19	CS-15192-SFX	Arrow 15 Fr. x 7.5" (19 cm) Radiopaque Polyurethane Catheter	19
23	CS-15232-SFX	Arrow 15 Fr. x 9" (23 cm) Radiopaque Polyurethane Catheter	23
27	CS-15272-SFX	Arrow 15 Fr. x 10.5" (27 cm) Radiopaque Polyurethane Catheter	27
31	CS-15312-SFX	Arrow 15 Fr. x 12" (31 cm) Radiopaque Polyurethane Catheter	31
	CS-15422-SFX	Arrow 15 Fr. x 16.5" (42cm) Radiopaque Polyurethane Catheter	42
	CS-15502-SFX	Arrow 15 Fr. x 19.7" (50cm) Radiopaque Polyurethane Catheter	50


Arrow® NextStep® Antegrade Chronic Hemodialysis Catheter Sets with SmartSeal® Hemostatic Sheath

PRODUCT NO.	DESCRIPTION	TIP-TO-CUFF INSERTION LENGTH (CM)
AC-15192-SFX	Arrow 15 Fr. x 7.5" (19 cm) Radiopaque Polyurethane Catheter	19
AC-15232-SFX	Arrow 15 Fr. x 9" (23 cm) Radiopaque Polyurethane Catheter	23
AC-15272-SFX	Arrow 15 Fr. x 10.5" (27 cm) Radiopaque Polyurethane Catheter	27
AC-15312-SFX	Arrow 15 Fr. x 12" (31 cm) Radiopaque Polyurethane Catheter	31
AC-15422-SFX	Arrow 15 Fr. x 16.5" (42cm) Radiopaque Polyurethane Catheter	42
AC-15502-SFX	Arrow 15 Fr. x 19.7" (50cm) Radiopaque Polyurethane Catheter	50

Arrow Simplicity® Micro-Introducer Kit

PRODUCT #	DESCRIPTION	QTY/BOX
MP-04450-AS	Color-Coded Micro-Introducer Set	5 ea.

Featured Kit Components

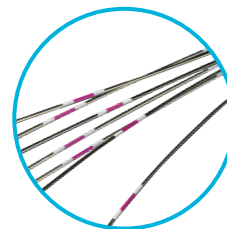

Arrow Simplicity® Color-Coded Micro-Introducer Kit

Unique system designed to simplify catheterization and facilitate optimal placement. Corresponds to color-coded packaging for Arrow® NextStep® Chronic Hemodialysis Catheters exposure to bloodborne pathogens.



Arrow SmartSeal® Hemostatic Peelable Dialysis Sheath

Designed to minimize risk of air embolism and clinicians' exposure to bloodborne pathogens.

Teleflex is a global provider of medical technologies designed to improve the health and quality of people's lives. We apply purpose-driven innovation – a relentless pursuit of identifying unmet clinical needs – to benefit patients and healthcare providers. Our portfolio is diverse, with solutions in the fields of vascular and interventional access, surgical, anesthesia, cardiac care, urology, emergency medicine and respiratory care. Teleflex employees worldwide are united in the understanding that what we do every day makes a difference. For more information, please visit teleflex.com.

Teleflex is the home of Arrow®, Deknatel®, Hudson RCI®, LMA®, Pilling®, Rüsch®, UroLift®, and Weck® – trusted brands united by a common sense of purpose.

Corporate Office

Phone +1 610 225 6800, 550 E. Swedesford Road, Suite 400, Wayne, PA 19087, USA

Regional Offices

United States: Phone +1 919 544 8000, Toll Free 866 246 6990, cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

Latin America: Phone +1 919 433 4999, la.cs@teleflex.com, 3015 Carrington Mill Boulevard, Morrisville, NC 27560, USA

International: Phone +353 (0)9 06 46 08 00, orders.intl@teleflex.com, Teleflex Medical Europe Ltd., IDA Business and Technology Park, Dublin Road, Athlone, Co Westmeath, Ireland

Does not contain natural rubber latex.
All Chronic Catheters sold individually.
Note: Products are not CE marked.

Caution: U.S. federal law limits this device to sale by or on order of a physician. Contents of unopened, undamaged package are sterile. Disposable. Refer to Instructions for Use (IFU) for current warnings, indications, contraindications, precautions and Instructions For Use.

Contact Teleflex for reference information on the above.

REFERENCES:

- 1 In vitro tests performed by an independent laboratory using 3:1 glycerine/saline at 37°C (Citech, Plymouth Meeting, PA).
- 2 Alcohol, alcohol-based solutions such as Chloraprep®; iodine-based solutions such as Povidone-iodine, mineral oil-based ointments, hydrogen peroxide or commercially available ExSept® are accepted for use with the Arrow® NextStep® catheter

Tegaderm is a registered trademark of 3M Company. Bactroban is a registered trademark of GlaxoSmithKline. Hibiclens is registered globally to one or more of the Mölnlycke Health Care Group of Companies. Chloraprep is a registered trademark of CareFusion. ExSept Plus is a registered trademark of Alcavis HDC, LLC or one of its affiliates.

Teleflex, the Teleflex logo, Arrow, Arrow Advancer, Arrow Simplicity, Deknatel, Hudson RCI, LMA, NextStep, Pilling, Rüsch, SharpsAway, SmartSeal, UroLift, and Weck, are trademarks or registered trademarks of Teleflex Incorporated or its affiliates, in the U.S. and/or other countries. Information in this material is not a substitute for the product Instructions for Use. Not all products may be available in all countries. Please contact your local representative. Revised: 05/2021.

© 2020 Teleflex Incorporated. All rights reserved. MC-006693 Rev 1